


Pressure Safe™

Neonatal High Flow Nasal Cannula



The New Measure
of Safety in Neonatal
High Flow Oxygen Delivery

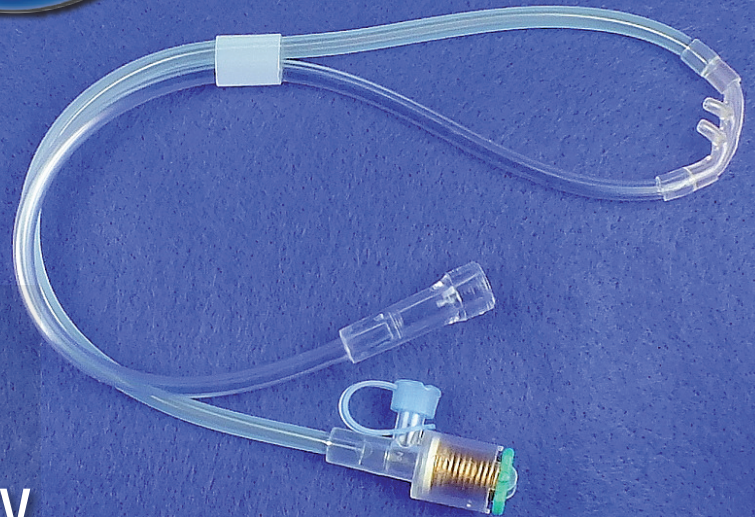


Westmed
When Your Care Is Critical

Pressure Safe™

High Flow Cannula for Neonatal Oxygen Therapy



The Pressure Safe High Flow Cannula is designed to monitor high pressure changes that may cause barotrauma conditions. A unique pressure limiting Pop-Off feature prevents airway pressure from exceeding 21 cm H₂O. The required pressure phases of high flow therapy are now clearly controlled and observed with Pressure Safe. In addition, the Comfort Soft Plus® patented elastomeric formulation minimizes cannula friction on the patient skin area. This produces a greater wear comfort and safeguards against epidermal stripping that can occur from traditional PVC cannula and/or adhesive securement tape use.

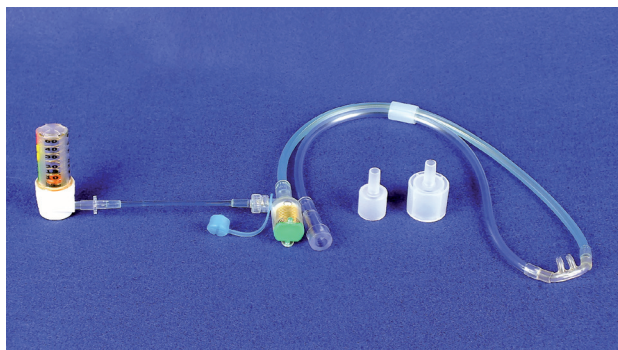

Compare These Advantages

- ✓ Pressure Pop-Off valve (will not exceed 21cm H₂O)
- ✓ Comfort Soft Plus® elastomeric material
- ✓ Sizes: Infant; Neonatal; and Large Infant
- ✓ Connection adaptors: 15mm OD; 15mm ID x 22mm OD
- ✓ Configured for multiple heated humidification systems
- ✓ Latex and DEHP Free

Item Code	Size/Description	Flow Rate	Quantity
0882	Neonatal Comfort Soft Plus® Nasal Cannula	Up to 10 LPM	10
0885	Infant Comfort Soft Plus® Nasal Cannula	Up to 10 LPM	10
0886	Infant w/manometer Comfort Soft Plus® Nasal Cannula	Up to 10 LPM	10
0887	Large Infant Comfort Soft Plus® Nasal Cannula	Up to 11 LPM	10

References:

Collins CL, Holberton JR, Barfield C, Davis PG. A randomized controlled trial to compare heated humidified high-flow nasal cannulae with nasal continuous positive airway pressure postextubation in premature infants. *J Pediatr*. 2013;162(5):949-954 e941.

Yoder BA, Stoddard RA, Li M, King J, Dirnberger DR, Abbasi S. Heated, humidified high-flow nasal cannula versus nasal CPAP for respiratory support in neonates. *Pediatrics*. 2013;131(5):e1482-1490.

Westmed
When Your Care Is Critical

Westmed, Inc.

5580 S. Nogales Highway
Tucson, AZ 85706
800-975-7987 | 520-294-7987
Fax: 520-294-6061
www.westmedinc.com

Distributed by:

CONNECT MEDICAL CORPORATION

118 Homestead Drive Suite 3

Phone: 1-866-416-NICU (6428) • Fax: 1-866-415-8962

www.connectmedical.biz • info@connectmedical.biz

U.S. Patent No. 8,333,200. US and International Patents Pending.
©2015 Westmed, Inc. All rights reserved. PN 75346 Rev .03



Visit www.westmedinc.com
for a complete listing of
Westmed products.